


# SternaLock<sup>®</sup> XP

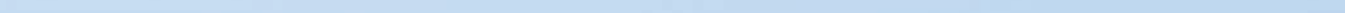
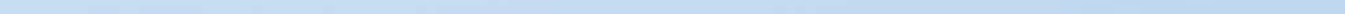
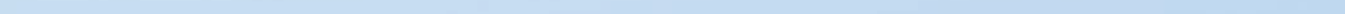
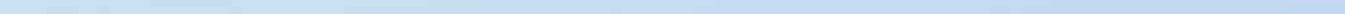
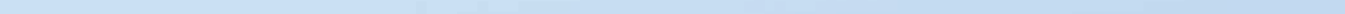
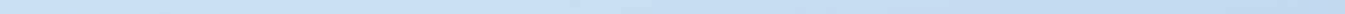
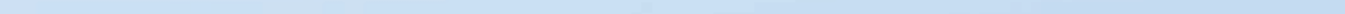
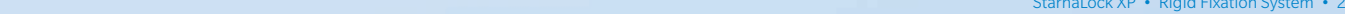
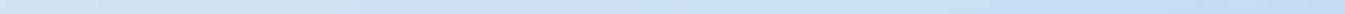
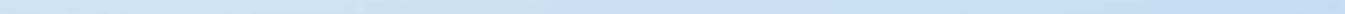
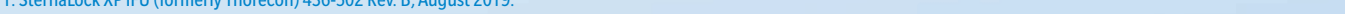
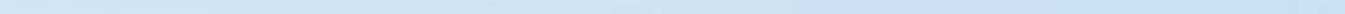


Rigid Fixation System

## SURGICAL TECHNIQUE

A medical device consisting of a central metallic plate with multiple circular holes and a central locking mechanism, connected to a thick, braided stainless steel cable. The plate is positioned in the center-right of the frame, with the cable looping around it. The background is a solid light blue.

# **SternaLock<sup>®</sup> XP**

## Rigid Fixation System



# SternaLock XP

Rigid Fixation System



# SternaLock XP Features

---

## **Ease of Use**

Disposable single use sterile system

## **Innovative Amorphous Cable**

Expands when tensioned to increase surface area<sup>1</sup>

## **Stronger and More Rigid Closure**

Designed to enhance stability and strength of traditional sternal closure technique<sup>2</sup>

## **Osteoporotic Indications**

Intended for use in patients with normal and/or poor bone quality<sup>2</sup>

---

1. Data on file - TR17063

2. SternaLock XP IFU (formerly Thorecon) 436-502 Rev. B, August 2019.

---

# **SternaLock XP**

## Surgical Technique

---

Surgical Technique shown is not applicable to the Tower Plate. Please reference the Mini-Thoracotomy Plate Surgical Technique (MKT940016).





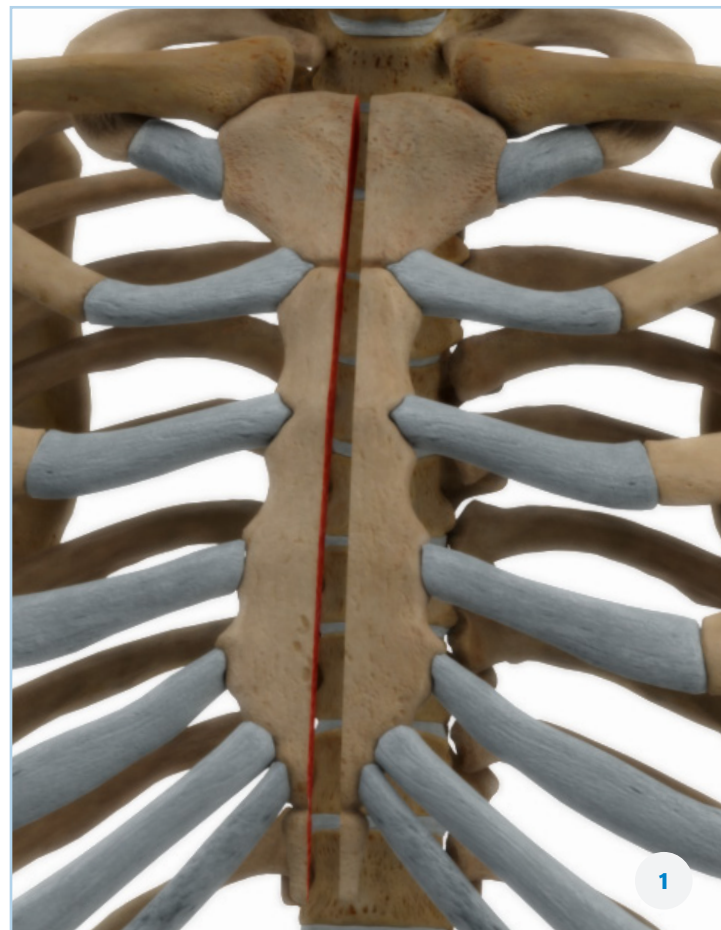
# 1. Exposure

Beginning medially, dissect all soft tissue from the surface of the sternum to allow for complete visualization; this may decrease the likelihood of an off-midline sternotomy. Optimally, the costal cartilage on both sides of the sternum is revealed. Perform a sternotomy (Figure 1).

## FOR SECONDARY CLOSURE OF THE STERNUM

Remove existing wires. Debride the involved sternal edges until they are free of devitalized tissues and hemostasis has been obtained. To allow for proper anatomical reduction and plate placement, bony calluses should also be removed from the midline and sternal surface (in addition to sternal soft tissue dissection in a revision surgery). A rib curette may be used to remove any nonviable cartilaginous rib tissue.

**Note:** A sternal bone specimen should be sent to pathology to assess for osteomyelitis and to microbiology for culturing. This may help the infectious disease consultant determine the appropriate antibiotic treatment if necessary.






## 2. Screw Verification

14mm screws are provided per procedural pack. Confirm 14mm screw lengths are appropriate by using the 14mm screw length indicator edge on the Torque Limiting Driver to measure the thickness of the sternal edge. (Figure 2). If leading edge of the torque driver protrudes past the posterior sternal edge, a 10mm screw should be used in this area of the sternum.

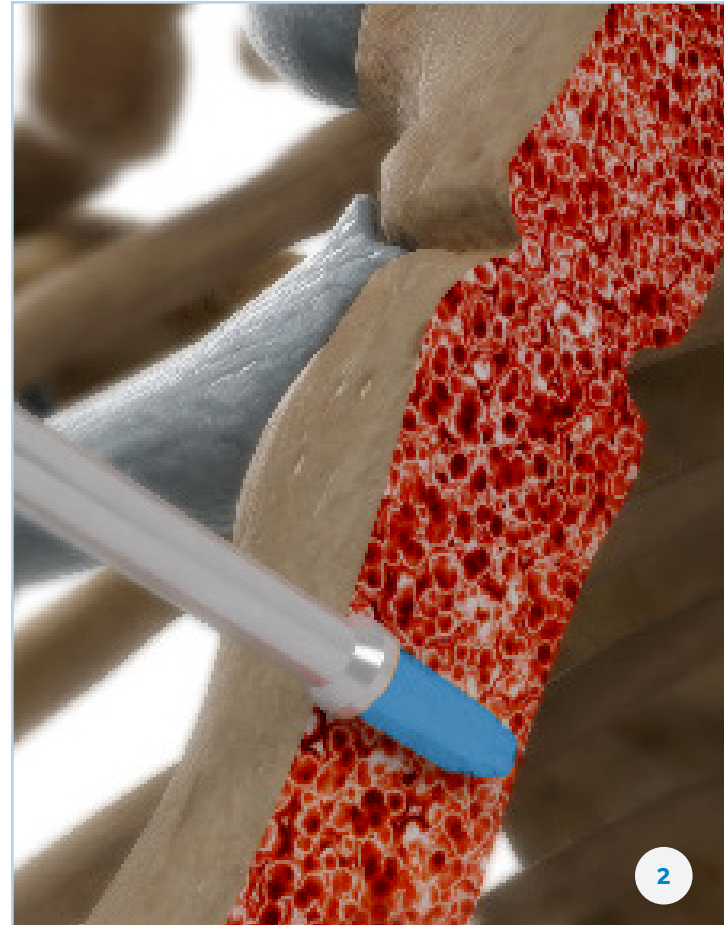Alternate screw options are available in separate packages, including 2.7mm (standard) and 3.0mm (rescue) diameters, in lengths of both 10mm & 14mm.





### 3. Cable Insertion

Once the plate is placed, route the cable around the sternum at the desired location. Cut the leader of the cable to remove the needle, leaving a small amount of leader to help guide the cable through the plate (Figure 3a). Thread the cable through the plate (Figure 3b).

**CAUTION:** Plates are not intended to be contoured or modified.

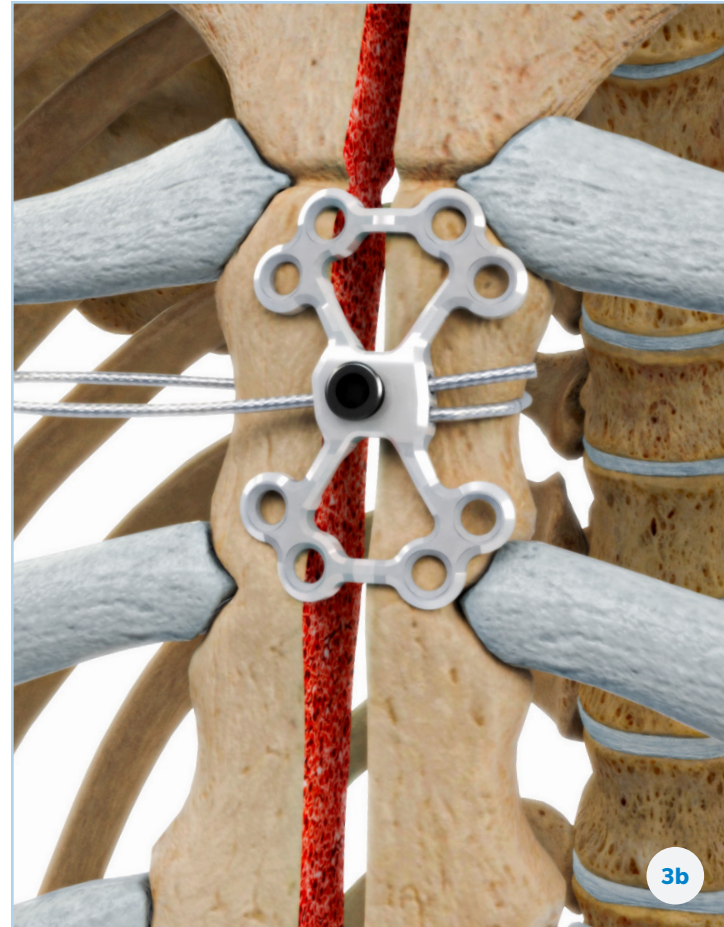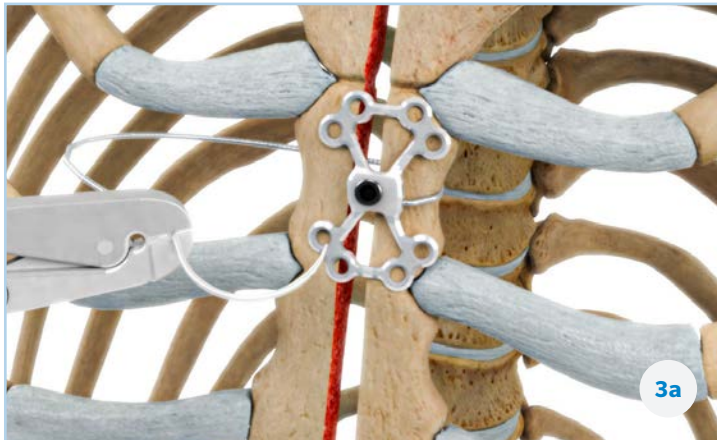**CAUTION:** When cutting the leader of the cable, there must be a small amount of leader remaining, otherwise the cable could fray and become difficult to pass through the plate and/or cable tensioner.



## 4. Sternal Approximation

---

The sternum can be reduced with the assistance of appropriate sized stainless steel sutures or stainless steel sternal cable at the xyphoid and manubrium. During reduction, observe the midline for protruding internal tissue and proper bony alignment. Mid-body sternal approximation can be completed using the included Cable Tensioner.

**Note:** Sternal approximation and stability must be determined by tactile assessment. Optimal tension is dependent on proper assessment of the patient's bone quality.

**CAUTION:** Putting dissimilar metals and alloys in contact with each other may be detrimental to the patient and/or function of the implant(s); therefore stainless steel cable or suture wire is recommended.

**WARNING:** Do not place any SternaLock XP plate over any other implants such as, but not limited to, plates, wires, screws or cables.

## 5. Cable Tensioning

To ensure a clear cable passage through the Cable Tensioner, rotate the handle clockwise until the line on the rear of the tensioner and on the handle are parallel. Pass the cable through the distal tip of the Cable Tensioner and advance the Cable Tensioner against the plate by pulling on the loose end of the cable to remove any slack. Obtain appropriate cable tension by rotating the handle clockwise (Figure 4).

**CAUTION:** Exercise caution while tensioning the cable. Over-tensioning could result in cable cutting into the bone, or cable breakage.

**Note:** Simultaneously tension the two cables that are most cephalad and caudal to better approximate the sternal constructs. This assures tension consistency between constructs and provides uniform tension across the sternum.



## 6. Bone Screw Placement

Using the Torque Limiting Driver's "stab-and-grab" feature, attach the bone screw to the Torque Limiting Driver (Figure 5a) and insert perpendicular through the desired screw hole on the plate. While applying downward force, rotate clockwise to drive the screw into the sternum (Figure 5b).





## 7. Screw Tightening

Fully seat the screw head into the plate. The screw head will thread into the plate causing the screw to lock to the plate (Figure 6).

**CAUTION:** Misalignment of the Torque Limiting Driver to the screw head may cause stripping of the screw head. Screws are not intended to be angled.

**Note:** Screw seating into the plate does not require the Torque Limiting Driver to reach full torque.

**Note:** Do not fully seat the screw into the plate when using a Power Driver. Fully seating the screws to the plate must be done manually using the Torque Limiting Driver.





## 8. Cerclage Locking

Insert the Torque Limiting Driver perpendicular into the Cerclage Set Screw. While applying downward force, rotate clockwise to drive the Cerclage Set Screw into the plate until the Torque Limiting Driver audibly clicks indicating desired torque (20in lbs.) is achieved and the Cerclage Set Screw is fully seated (Figure 7).



## 9. Cable Cutting and Tension Removal

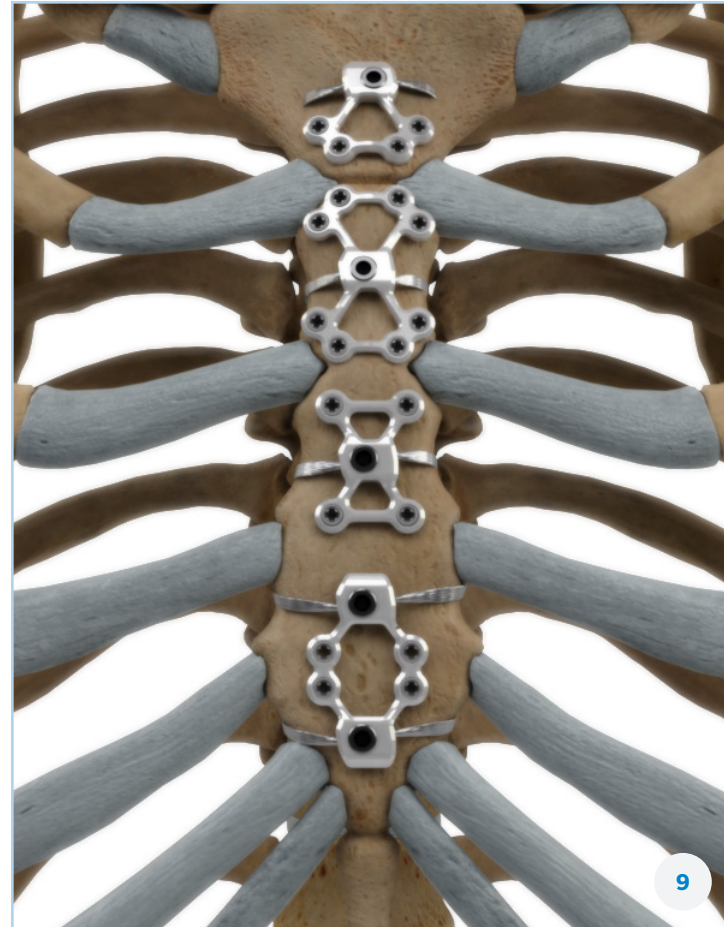
Fully insert the Torque Limiting Driver into the hollow tube on the distal tip of the Cable Tensioner (Figure 8a). Gently lever the handle forward until the Cable Tensioner releases from the plate indicating the cable is fully cut flush with the plate (Figure 8b). Once the cable has been cut, remove the Cable Tensioner, Torque Limiting Driver and discard along with remaining cable (Figure 8c).

**CAUTION:** The cable must only be cut once, flush to the plate, and prior to removing the Cable Tensioner from the device.



## 10. Complete Construct

It is important to be aware of the overall stability of the closure and use as many plates as necessary to achieve adequate fixation based on a surgeon's assessment for each patient (Figure 9).



# 11. Closure

---

After the plate/screw constructs, and all cable or wire (if applicable) are installed in final position, wound closure is performed. Close skin and subcutaneous tissues in standard closing procedure.

**Note:** Do not pull or lift the patient by the arms. It is also not recommended to raise the arms higher than 90° at shoulder level.

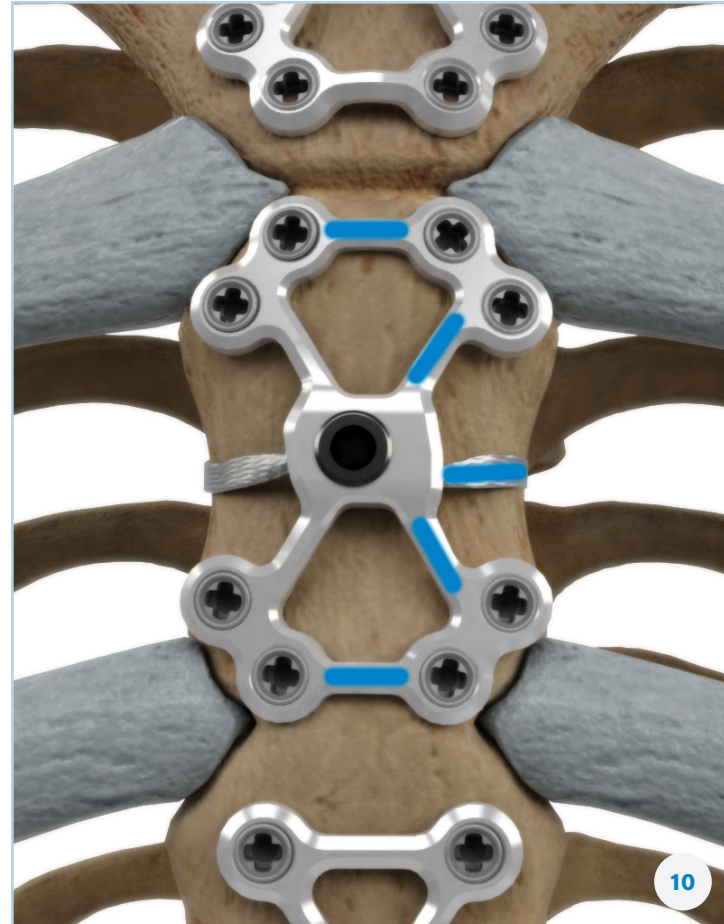
## 12. Emergent Re-entry/Removal

### Extraction (if required)

To remove the bone screws, insert the Torque Driver into the head of the screw and rotate counter-clockwise while exerting downward pressure. Use a cable cutting instrument to cut the cable from the sternum.

### Emergent Re-Entry (if required)

First, use a cable cutting instrument to cut the cable. Next, cut the plate through the main support locations (Figure 10; 8 Hole X-Plate shown as example). If needed, follow extraction instructions above.
















## SternaLock XP • Implants

<p><b>94-1200-08:</b> X-Plate Kit</p> <p>8x 14mm Screws, 1x Cable, 1x Tensioner, 1x Driver</p>	
<p><b>94-1300-04:</b> Double Cable Plate Kit</p> <p>4x 14mm Screws, 2x Cable, 2x Tensioner, 1x Driver</p>	
<p><b>94-1300-04-S:</b> Double Cable Plate Kit (Sharp)</p> <p>4x14mm Screws, 2x Cable, 2x Tensioner, 1x Driver</p>	
<p><b>94-1400-04:</b> Box Plate Kit</p> <p>4x 14mm Screws, 1x Cable, 1x Tensioner, 1x Driver</p>	
<p><b>94-1500-04-S :</b> Auxiliary Plate Kit</p> <p>4 x 14mm Screws, 1x Cable,1x Tensioner, 1x Driver</p>	
<p><b>94-1600-05:</b> Tower Plate Kit</p> <p>2x 10mm Screws, 3x 14mm Screws, 1x Cable, 1x Tensioner, 1x Driver</p>	
<p><b>94-1700-10:</b> Ladder Plate Kit</p> <p>10x 14mm Screws, 2x Cable, 2x Tensioner, 1x Driver</p>	


Items not to scale

## SternaLock XP • Optional Instrumentation

<p><b>94-Driver-S:</b> Driver</p> <p>Ordering Qty: 1</p>		<p><b>70-2001:</b> ZDrive™ High Torque Driver &amp; Battery</p> <p>Ordering Qty: 1</p>	
<p><b>94-Tensioner-S:</b> Cable Tensioner</p> <p>Ordering Qty: 1</p>			
<p><b>94-Platecutter:</b> Plate and Cable Cutter</p> <p>Ordering Qty: 1</p>		<p><b>73-1195:</b> ZDrive™ Blade, SternaLock XP</p> <p>Ordering Qty: 1</p>	

Items not to scale

## SternaLock XP • Additional Sterile Screws

Part Number	Description	Content	
94-27-14-8S	Bone Screw 2.7 x 14mm	8 Pack	
94-27-10-8S	Bone Screw, 2.7 x 10mm	8 Pack	
94-27-14-6S	Bone Screw 2.7 x 14mm	6 Pack	
94-27-10-6S	Bone Screw, 2.7 x 10mm	6 Pack	
94-27-14-4S	Bone Screw 2.7 x 14mm	4 Pack	
94-27-10-4S	Bone Screw, 2.7 x 10mm	4 Pack	
94-30-14-2S	Bone Screw, 3.0 x 14mm	2 Pack	
94-30-10-2S	Bone Screw, 3.0 x 10mm	2 Pack	

Items not to scale



**For more information on the SternaLock® XP Rigid Fixation System and other solutions, please contact us at:**

**BIOMET MICROFIXATION GLOBAL HEADQUARTERS**

1520 Tradeport Drive • Jacksonville, FL 32218-2480 • Tel 904.741.4400 • Toll-Free 800.874.7711  
Fax 904.741.4500 • Order Fax 904.741.3059 • [zimmerbiomet.com](http://zimmerbiomet.com) • [zbthoracic.com](http://zbthoracic.com)

SternaLock XP Rigid Fixation System is manufactured by

**A&E MEDICAL CORPORATION**

5206 Asbury Road, Farmingdale, NJ 07727 • Tel: 904.741.4400 • [aemedical.com](http://aemedical.com)

All content herein is protected by copyright, trademarks and other intellectual property rights owned by or licensed to Zimmer Biomet, or its affiliates unless otherwise indicated, and must not be redistributed, duplicated or disclosed, in whole or in part, without the express written consent of Zimmer Biomet. This material is intended for healthcare professionals and the Zimmer Biomet sales force. Distribution to any other recipient is prohibited. Zimmer Biomet does not practice medicine. Because this information does not purport to constitute any diagnostic or therapeutic statement with regard to any individual medical case, each patient must be examined and advised individually, and this information does not replace the need for such examination and/or advice in whole or in part. Each physician should exercise his or her own independent judgment in the diagnosis and treatment of an individual patient. Check for country product clearances and reference product-specific instructions for use. For complete product information, including indications, contraindications, warnings, precautions and potential adverse effects, see the package insert, or visit the manufacturers website.